

GTR-ID

Industrial Tube Rotator



GTR-ID Industrial Tube Rotator shown with included 50mL tube holder disk

The Globe Scientific GTR-ID Industrial Tube Rotator has a rugged, workhorse quality and a value proposition that simply can't be beat!

The GTR-ID Tube Rotator efficiently mixes your biological samples and solutions using the included 50ml tube holder disk or the optional 1.5mL/2.0mL and 15mL disks. An industrial-style metal housing with tilt adjustment ensures a lifetime of worry-free operation. Consistent tilt adjustment is performed by referring to the graduated angle scale, which is printed directly on the unit housing.

The unit housing, tube holder disks and tube clamps and are metal, alleviating the breaking of plastic holders and clamps seen in models from other manufacturers. Optional disk stacking rods allow up to two tube disk holders to be used simultaneously for increased productivity.

Features:

- Heavy-duty metal chassis and tube holders/disk accessories
- Digital adjustment of speed and time
- LCD display
- Speed range from 10 to 70rpm
- Full adjustment of mixing action from end over end or rolling
- Adjustable mixing angle from 0° to 90°
- A wide variety of accessories are available
- Cold room friendly
- Can stack up to two tube holders/disks for additional capacity


Industrial Tube Rotator



Digital display of the GTR-ID can be accurately set using the parameter adjustment knob



Included tube holder disk accommodates (8) 50mL tubes – additional tube holder disks are available

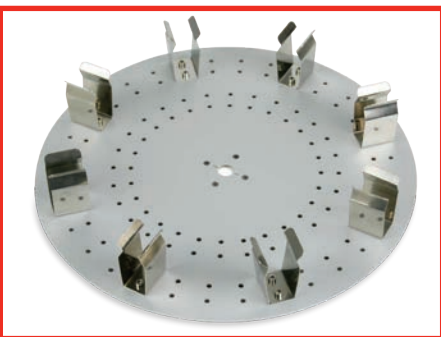

Robust stainless steel tube holder disks for years of trouble free use

GTR-ID Specifications

Tilt Angle	0° - 90°
Speed Range	10 - 70rpm
Speed Display	LCD
Timer Display	LCD
Timer	1-1199min
Operation Mode	Continuous or timed
Motor Type	DC motor
Voltage	100-240V, 50/60Hz
Power	40W
Dimensions (WxDxH)	280 x 210 x 300mm (11.0 x 8.3 x 11.8in)
Weight	4.5kg (10lbs)
Operating Temperature	5° - 40° C
Permissible Relative Humidity	80% RH (non-condensing)
Protection Class	IP21
Warranty	2 Years

ORDERING INFORMATION

Item #	Description	Unit
GTR-ID	Industrial Tube Rotator, Digital, Variable Speed and Angle, 120-240v, 50/60Hz, (Includes: 8-Place Disk for 50mL Tubes)	Each
GTR-ID-1.5	Tube Holder Disk, 60-Place, for 1.5mL Microcentrifuge Tubes	Each
GTR-ID-15	Tube Holder Disk, 16-Place, for 15mL Centrifuge Tubes	Each
GTR-ID-50	Tube Holder Disk, 8-Place, for 50mL Centrifuge Tubes	Each
GTR-ID-RODS	Tube Holder Disk Stacking Rods (allows use of 2 GTR-ID Tube Holder Disks at once)	4

www.globescientific.com

Customer Service, Sales and Technical Support:

Phone 1 (800) 394-4562 • 1 (201) 599-1400

Fax 1 (201) 599-1406

E-mail mail@globescientific.com

globe[®]
SCIENTIFIC INC.